



Life Changing Spinal Innovation



HORIZON

The Right Size for Every Body



Offering Effective Support and Pain Relief

ADJUSTABLE Horizon Braces



Horizon™ 456 TLSO
CODE L0456 APPROVED

Horizon™ 637 LSO
CODE L0637 APPROVED

Horizon™ 631 LSO
CODE L0631 APPROVED

Horizon™ 627 Lumbar
CODE L0627 APPROVED

The Horizon™ line of fully adjustable, inelastic spinal braces offers a clinically proven treatment option for patients experiencing low back pain. Combining a therapeutic level of support and direct independent compression, Horizon braces provide a non-invasive, non-narcotic solution to relieve pain and enable patients to regain activities of daily living.

The Right Size for Every Body

24-70_{IN}

Simple pull-through sizing allows one brace to fit waists ranging from 24 to 50 inches

and up to 70 inches with the extension panels. Horizon's ergonomic and low-profile design may be comfortably worn under clothing.

Minimum Inventory, Maximum Cost Savings

One size adjustable means the correct size is always on hand, minimizing inventory requirements and significantly reducing costs.



Inelastic Braces Proven Effective in Controlling Low Back Pain

Research shows that inelastic LSO's are significantly more effective at improving trunk stability than elastic braces,^{1,2,3} a key factor in relieving pain. Increased trunk stability reduces the load on overtaxed muscles,⁴ which can dramatically reduce low back pain. Aspen braces are inelastic by design for maximum pain relief.

1. Morrisette DC, Cholewicki J, Logan S, Seif G, McGowan S. A randomized clinical trial comparing extensible and inextensible lumbosacral orthoses and standard care alone in the management of lower back pain. *Spine*. 2014.
2. Cholewicki J, Lee AS, Reeves NP, Morrisette DC. Comparison of trunk stiffness provided by different design characteristics of lumbosacral orthoses. *Clinical Biomechanics*. 2010.
3. Cholewicki J, Reeves NP, Everding VQ, Morrisette DC. Lumbosacral orthoses reduce trunk muscle activity in a postural control task. *Journal of Biomechanics*. 2007.
4. van Dieën JH, Selen LPJ, Cholewicki J. Trunk muscle activation in low-back pain patients, an analysis of the literature. *Journal of Electromyography and Kinesiology*. 2003.



Mold • Bend • Trim Easy Access for Custom Modification

All structural components are easily removed and can be heat molded, bent or trimmed to further customize the brace to accommodate individual patient conditions.



One Size Adjustable

Pull-through sizing mechanism allows for quick and easy adjustment, without the use of tools.

Integrated Anterior Panels

Anterior panels at each end of the brace overlap to guarantee effective midline support. Aspen's innovative FlexTab™ Technology disperses pressure and naturally conforms to each patient's anatomy for increased comfort.

Horizon 637 Lateral Panels

Provides increased lateral control for patients in need of additional support. Panels can be removed to support patient recovery.



*Horizon 637 LSO shown in smallest configuration

Reliable Structural Integrity

Tri-band technology system and integrated vertical stays provide stability and structural integrity for effective compression, while allowing for flexibility where needed.


Stability and Comfort

Breathable mesh spacer fabric for increased comfort.



Custom Support

When tightening, firm but flexible back panel naturally conforms to the patient's lordosis to provide true structural support.

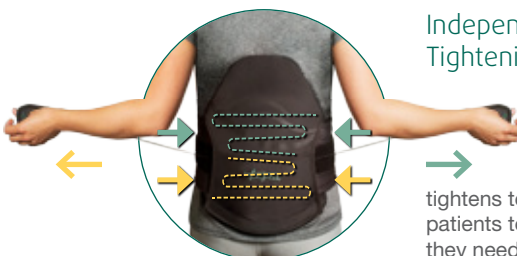

Independent Tightening Mechanism

Aspen's patented 4:1 mechanically advantaged SlickTrack™ tightening system independently tightens top and bottom allowing patients to localize compression where they need it for maximum pain relief.



A Height for Every Body

Telescoping shoulder straps on the Horizon 456 provide 13 inches of height adjustment to effectively fit even the tallest patients. This, along with circumferential adjustability, enables the right fit for every body.



Ergonomically Designed Molded Pull Tabs

Easy to locate by feel and designed for a comfortable pull. Bonded and stitched materials ensure their strength and durability.

The Right Brace for Every Condition

Regardless of which configuration a patient requires initially, all Horizon braces can be stepped down to the Horizon PRO (Pain Relief Orthosis), the lowest-profile option to help manage recurring back pain as needed.



Horizon PRO



The Aspen Advantage

Aspen Medical Products builds the finest spinal braces on the market. With an unparalleled commitment to Research and Development, Aspen has brought to market numerous unique features, marked by the Aspen Advantage logo. These features yield real benefits for practitioners and patients. In addition, Aspen braces are backed by more than 30 years of clinical research that substantiates the safety,

efficacy, superior comfort and improved patient outcomes that these products provide. Quality products you can count on with features that make a positive difference. That's the Aspen Advantage.




The Aspen Advantage Symbol

An innovative feature available only from Aspen

Giving Experts Modification Options

Aspen Spinal Braces provide unique opportunities to access structural components to remove, trim, bend, mold, assemble and provide individualized customization options. This enables experts to significantly modify each brace to treat specific patient indications, aid with the treatment of different diagnoses and adapt to different bracing treatment philosophies or physician direction. Made of high quality materials such as nylon, high-density and low-density polyethylene plastics and non-ferrous metals, these durable and resilient components not only provide the best support, they allow for significant modification without breaking, buckling or cracking that can occur with lower quality materials. Trimming, bending, and molding are only a few ways in which Aspen braces can be significantly modified to further improve patient outcomes.

For more information on how Aspen braces can be significantly modified, or for documentation templates, contact your local Aspen Sales Representative or visit [aspenmp.com/product-modifications](https://www.aspenmp.com/product-modifications).

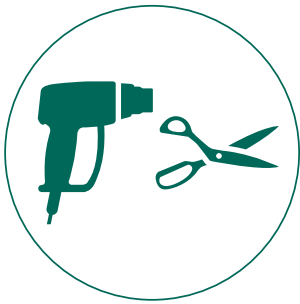

Sizing Information

Waist Circumference	
24 —————> 50	—————> 70
24-50 in (61-127 cm)	Up to 70 in (178 cm)
Horizon One Size Adjustable Horizon 456 TLSO · Horizon 637 LS0 · Horizon 631 LS0 · Horizon 627 Lumbar	Horizon Extension Panels Adds 20 in (51 cm)

Product Information

Horizon	Horizon 456 TLSO	993740
	Horizon 637 LS0	993730
	Horizon 631 LS0	993720
	Horizon 627 Lumbar	993710
Accessories	Horizon Extension Panels	993745
	Therapy Pack	990005